


Model Information & Spec



Control Box

Model Name		SWCLD-EP-LN1		SWCLD-EP-LN2		SWCLD-EP-SN1		SWCLD-EP-SN2		SWCLD-EP-SN3	
1	Epidural Catheter	SWCLD-EC-05		SWCLD-EC-05		SWCLD-EC-02		SWCLD-EC-02		SWCLD-EC-06	
		Length (mm)	OD(Ø)	Length (mm)	OD(Ø)	Length (mm)	OD(Ø)	Length (mm)	OD(Ø)	Length (mm)	OD(Ø)
		200	4.1	200	4.1	120	4.1	120	4.1	100	4.1
2	Shaver	SWCLD-ET-07		SWCLD-ET-06		SWCLD-ET-03		SWCLD-ET-05		SWCLD-ET-08	
											
		Length (mm)	OD(Ø)	Length (mm)	OD(Ø)	Length (mm)	OD(Ø)	Length (mm)	OD(Ø)	Length (mm)	OD(Ø)
		320	3.5	320	2.5	235	3.5	235	2.5	215	3.5
3	Dilator	SWCLD-DI-04		SWCLD-DI-04		SWCLD-DI-03		SWCLD-DI-03		SWCLD-DI-05	
4	Surgical Probe	SWCLD-HK-03		SWCLD-HK-03		SWCLD-HK-01		SWCLD-HK-01		SWCLD-HK-04	
5	Trocar	SWCLD-ND-03		SWCLD-ND-03		SWCLD-ND-02		SWCLD-ND-02		SWCLD-ND-04	
6	Control Box	SWCLD-CB-01									

Epidural Catheter

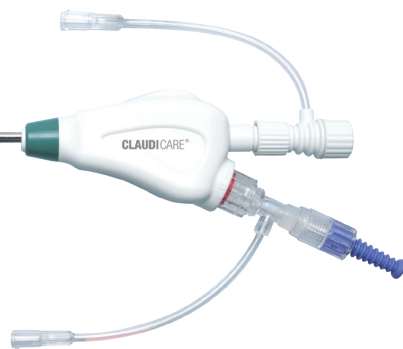

Shaver




CLAUDICARE®

Percutaneous
Neurolysis and Foraminoplasty



Trocar

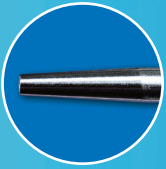

Dilator




Surgical Probe

Percutaneous Neurolysis and Foraminoplasty

Percutaneous Neurolysis and Foraminoplasty is the minimally invasive surgery for patients who have claudication or backpain resulting from foraminal stenosis. This surgery is used to widen narrowed foramen intervertebrale by approaching the immediate proximity of diseased part via the transforaminal route and removing the lesion (ligament, disc, soft tissue, etc.) through the minor incision (3mm) under local anesthesia.

Indication


- Lumbar foraminal stenosis
- Lumbar central stenosis
- Chronic backpain
- Patients refusing operation or inoperable patients
- Effective diagnostic root block but no longterm relief

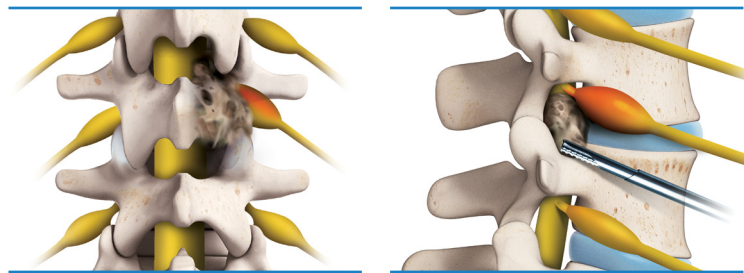

Benefits




Safe : Blunt micro bit & protection shield
Micro: 3mm external diameter bit
Easy : about 20 minute surgical time & fast recovery
 All-In-One Kit including a battery pack

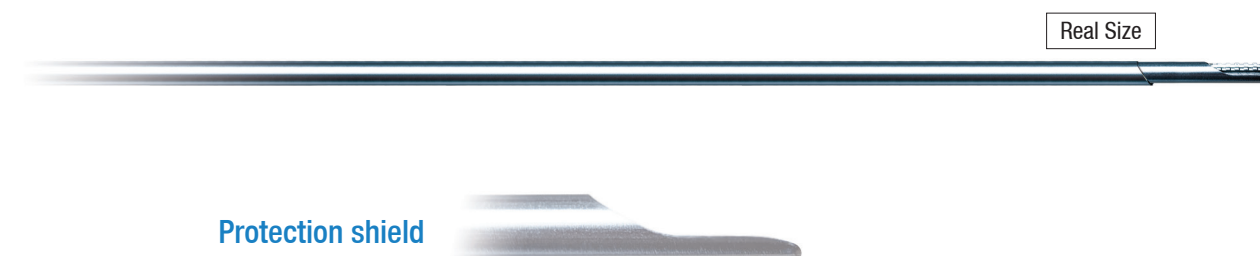


Foraminoplasty with claudicare



I Guide wire to Kambin's triangle

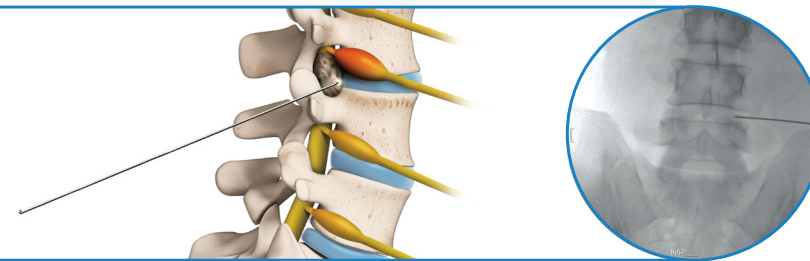

II Dilator & Working channel

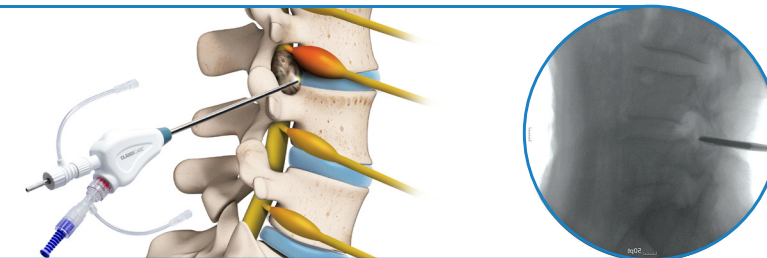

III Foraminoplasty with claudicare

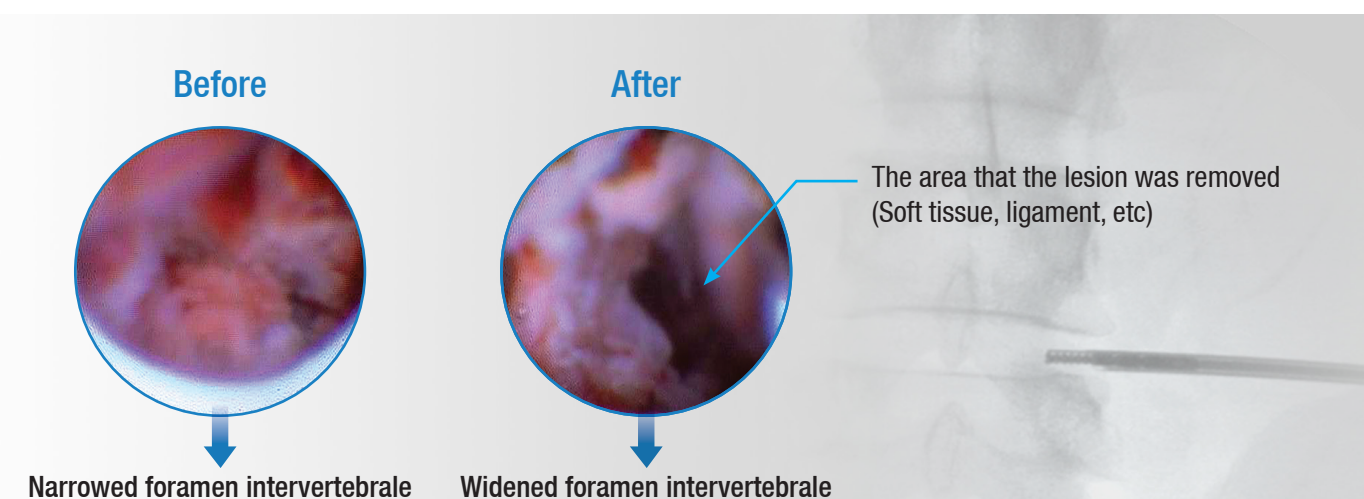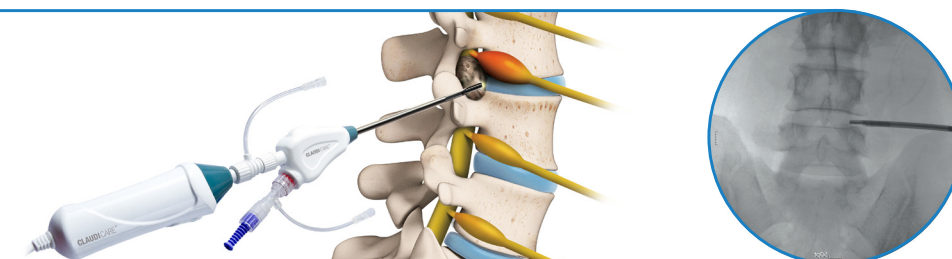

Model Information & Spec



Model Name		SWCLD-EP-02		SWCLD-EP-04		SWCLD-EP-S		SWCLD-EP-L	
1	Epidural Catheter	SWCLD-EC-02		SWCLD-EC-04		SWCLD-EC-02		SWCLE-EC-05	
		Length (mm)	OD(Ø)	Length (mm)	OD(Ø)	Length (mm)	OD(Ø)	Length (mm)	OD(Ø)
		120	4.1	120	3.4	120	4.1	200	4.1
2	Shaver	SWCLD-ET-01		SWCLD-ET-02		SWCLD-ET-03		SWCLD-ET-04	
									
		Length (mm)	OD(Ø)	Length (mm)	OD(Ø)	Length (mm)	OD(Ø)	Length (mm)	OD(Ø)
		260	3.5	260	2.5	235	3.5	320	3.5
3	Dilator	SWCLD-DI-01		SWCLD-DI-02		SWCLD-DI-03		SWCLD-DI-04	
4	Surgical Probe	SWCLD-HK-01		SWCLD-HK-02		SWCLD-HK-01		SWCLD-HK-03	
5	Trocar	SWCLD-ND-01		SWCLD-ND-01		SWCLD-ND-02		SWCLD-ND-03	
6	Battery Pack	SWCLD-BP-01							